



Dental Floss Series



DF201



DF201M(mint)



DF202



DF202M(mint)



DF203-W(white)



DF203-G(green)



DF203-O(orange)



DF203-B(blue)



DF204-P(pink)



DF204-G(green)



DF204-B(blue)




DF204-V(violet)




Printed card package

Usage Direction:

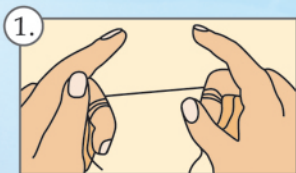

1. Break off 45cm of floss, wind the floss around the middle fingers of each hand.

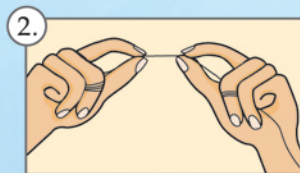

2. Tighten the floss with thumbs and forefingers, leaving 5cm of floss between them.



3. Gently guide the floss between teeth and run carefully up and down the side of each tooth to clean off plaque and debris. Repeat the above process for all teeth. Use a clean section of floss each time.